


Digital Picture

The digital picture will be overlaid on the CT scan data in an area as pictured.

- Picture taken face forward in scan position
- Lips in a closed and relaxed position
- Keep identified area clear of hair or foreign objects for best results

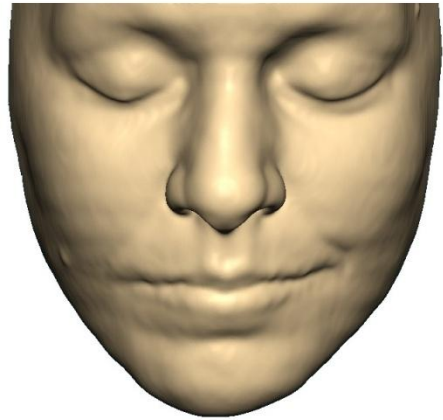

Facial area for soft tissue prediction

Make a CD or disc with the uncompressed DICOM image data.

Label the CD or disc with the following:

- Patient Identification
- Physician Name
- Facility Name
- CT Technologist Name
- Technologist Phone Number

You may upload the uncompressed DICOM data and digital picture to our secure ftp site or send them to the address below. You may also arrange an alternate secure data transfer method of your choice.

ProtoMED, Inc.

1329 West 121st Avenue

Westminster, CO 8023

Please contact us for instructions on uploading CT data or with any questions



Phone: 303-466-5610

Fax: 303-255-5576